

SonoScape

S22



Bring Ultrasound Excellence to Anywhere



SonoScape's S22 is a fully-featured ultrasound system in a slim and elegant package. It combines mobility with utility to fit in almost any clinical situation.

A host of proven and renowned technologies make such a difference, the clarity of the image improves accuracy and precision. The S22 is aimed at helping the patients and doctor feel more comfortable physically and emotionally during medical examination

S22

Brilliant Performance



With high definition resolution and unmatched penetration, S22 helps you improve the clinical confidence by affording exceptional image quality especially in detecting abnormalities or small lesion. Furthermore, intuitive protocols will simplify image acquisition process and enhance the overall image quality. They empower you to enhance spatial and contrast resolution and sharpen border definition for outstanding image clarity. Focus on image quality, S22 fulfills the utmost value during diagnosis

Phase Inversion Harmonic Imaging (PIH)

The harmonic signals are fully preserved without degradation of the acoustic information, which makes it possible for S22 to image high-level details and improved contrast resolution by reducing noise and clutter in the visualizing of subtle lesions, small parts, vascular and so on.

Special Compound imaging (SCI)

Spatial Compound Imaging utilizes several lines of sight for optimal contrast resolution, speckle reduction and border detection, with which S22 is ideal for superficial and abdominal imaging with better clarity and improved continuity of structures

Speckle Reduction Imaging (μ -scan)

The second-generation μ -scan technology greatly improves the visibility of organs and lesions with improved, high-definition contrast resolution that suppresses speckle artifact while maintaining real tissue architecture.

S22 boast its slender contour with cohesive and intuitive solution package which enable doctors to focus more on the patient, and much less on operation steps. Optimized file management, innovative console layout and SonoScape's uniquely renowned transducer technologies address every aspect of comfort for both doctors and patients.



Brand-new and powerful platform

Evaluating patients in today's demand in clinical environment, aim to reduced time and bring doctor more confidence, compact with the latest technology and optimized work flow, Sonoscape release the latest platform, can supply with abundant and dedicated clinical tools:

- Full range of transducers
- Customizable menu
- Enhanced work flow

Slim and ultra-compact design


Control Panel
18.5'inch flat LED Monitor
The latest backlight monitor, use the highest standard, try to provide the most contrast image and better color expression performance.

- Simple user interface


Light mobility, easy control

- Bar scanner
- Remote Control
- Wireless networking technology


Superior Imaging Quality

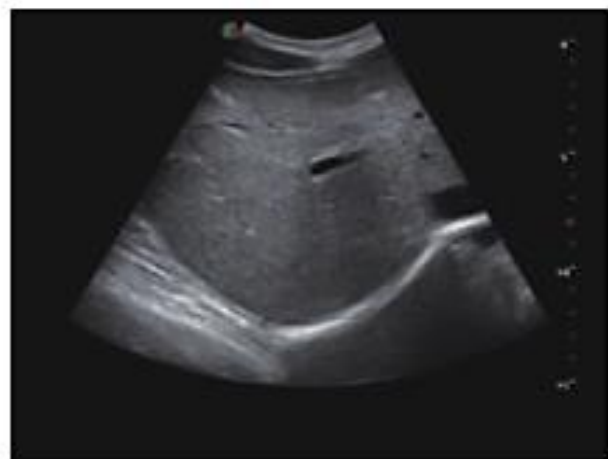


2D Liver

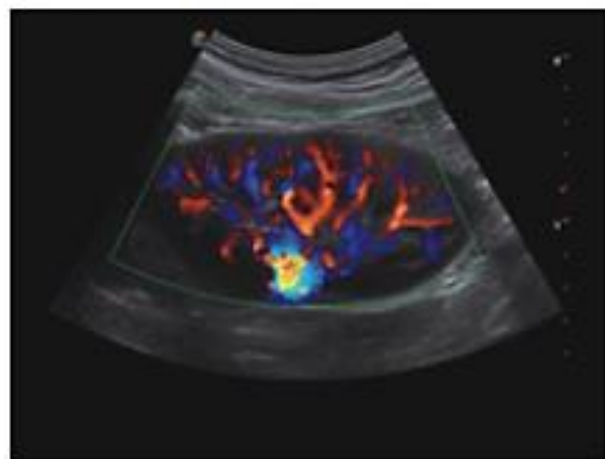


Kidney Flow

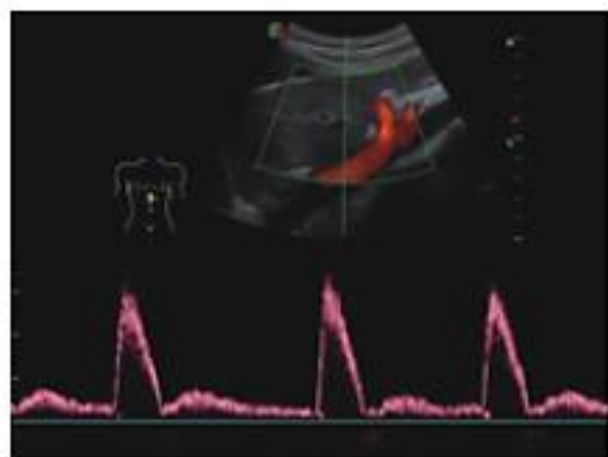

Abdominal Aorta PW

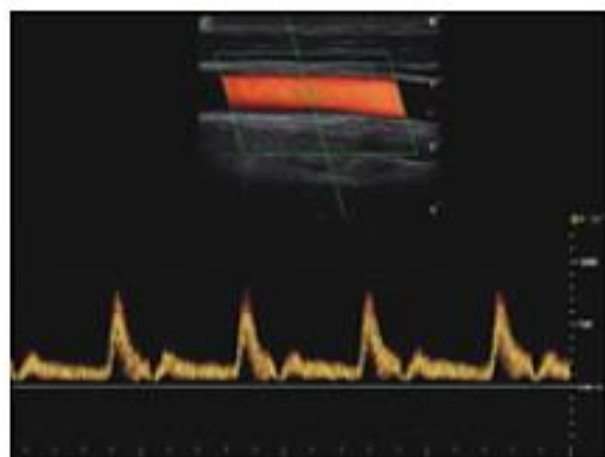

Common Carotid Artery PW

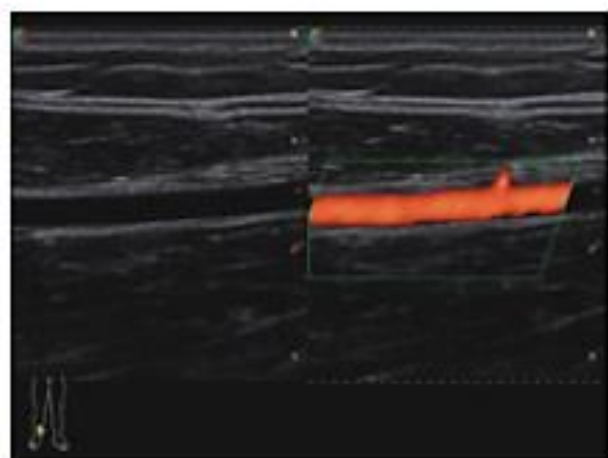

Dual Dynamic Display of Posterior Tibial Artery



Fetus